


SAVE TIME & MONEY!

The first vacuum separation tank that automatically drains when needed without the added cost of a pump.

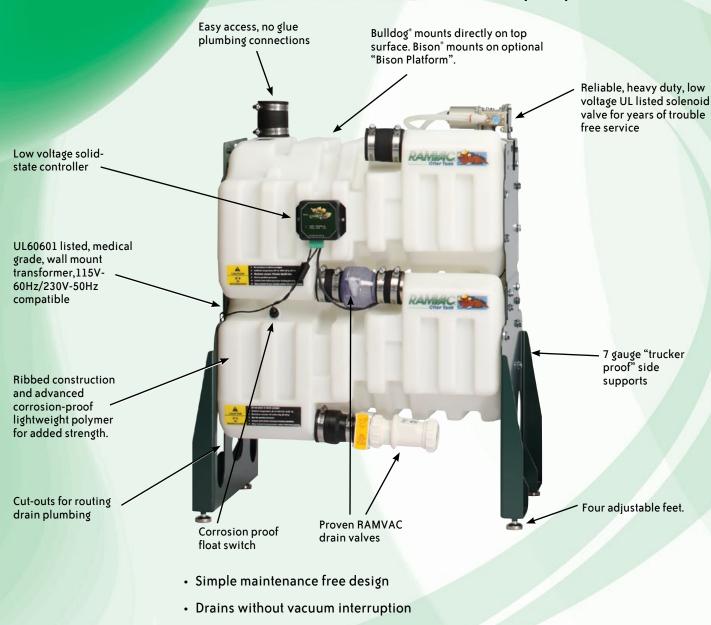

· Significantly reduces the size and number of separating tanks

Completely automated to provide worry free operation

· Intelligent controller records cycle history

· Fits neatly next to walls and into corners





What does the InfiniTank do and why do you need one? Most dry dental vacuum systems direct liquids and solids generated in the treatment room to a collection tank where they are stored until the vacuum unit is shut off. At that time, they are allowed to flow to a sanitary drain. In a properly designed system, this process happens once a day (usually at the end of each work day). If the separating tank has been oversized, there is wasted space and money. However, if the separating tank has been sized too small, your problems are even greater. You will have the nuisance of cycling your vacuum pump on and off during the day, resulting in irregular service and wasted time. You also run a constant risk of damaging your vacuum pump if the tank accidentally overfills.

The InfiniTank is the first vacuum separation tank that automatically drains without the use of expensive, maintenance prone pumps. The InfiniTank works in harmony with the natural forces of gravity and air pressure to create "infinite" tank volume in an easy to install, compact design. In short, there's no more cycling your vacuum pump during work hours and no more risk of flooding your vacuum pump because your separating tank is too small. Also, the InfiniTank can eliminate the need for mulitple separation tanks and as with other RAMVAC tank designs, you can mount the vacuum unit to the top of the InfiniTank to create an even more compact installation.

InfiniTank Part Number 575080 or 575085 (50Hz)

Installation Options

Standard InfiniTank

Place Bison in any convenient location. Tank does not have to be close by.

Dimensions: Width 29" (740 mm) Depth 19" (540 mm) Height 36" (920 mm)



InfiniTank with height extender Height Extender (575212)

Assembled Dimensions:

Width 29" Depth 19" Height 48"



InfiniTank with Bison® platform and Bison

Dimensions (w/ Bison):

Width 34" Depth 21" Height 70"



InfiniTank with Bison placed "side by side"

Dimensions (w/ Bison):

Width 55" Depth 22" Height 36"





InfiniTank with additional 50-gallon tank for high flow - high volume applications

Dimensions (w/ Bison):

Width 55" Depth 22" Height 62"





Patent Pending

Description Item# InfiniTank 575080

575085 InfiniTank 50 Hz



